



Mouse anti Tumor-Associated Glycoprotein (TAG-72) Monoclonal Antibody

Alternative Name(s): Tumor-Associated Glycoprotein (TAG-72)

Order Information

- **Description:** Tumor-Associated Glycoprotein (TAG-72)
- **Catalogue:** 603-880
- **Lot:** See label
- **Size:** 100ug/200ul
- **Host:** Mouse
- **Clone:** B72.3
- **Application:** IHC(P), WB
- **Reactivity:** Hu

ANTIGEN PREPARATION

A membrane extracts of human hepatic carcinoma

BACKGROUND

Tumor-Associated Glycoprotein (TAG-72) is a high-molecular-weight mucin, protein/sugar complex found on the surface of many cancer cells, including breast, colon, and pancreatic cells. It is a target for several anti-cancer drugs and a tumor marker. In studies TAG-72 showed a higher positivity rate for gastric cancer than CEA, CA 19-9, and AFP.

PURIFICATION

The mouse IgG is purified by Protein A-Affinity Chromatography according to Isotyping

FORMULATION

This affinity purified antibody is supplied in sterile Phosphatebuffered saline (pH7.2) containing antibody stabilizer

SPECIFICITY

This antibody recognizes human Tumor-Associated Glycoprotein (TAG-72) protein. The other species are not tested.

STORAGE

The antibodies are stable for 24 months from date of receipt when stored at -20°C to -70°C. The antibodies can be stored at 2°C-8°C for three months without detectable loss of activity. Avoid repeated freezing-thawing cycles.

APPLICATIONS/SUGGESTED WORKING DILUTIONS*

- Western Blot: 0.1-1 µg/ml
- ELISA: 0.01-0.1 µg/ml
- Immunoprecipitation: 2-5 µg/ml
- IHC: 2-10 µg/ml
- Flow cytometry: Not tested
- Molecular Weight: 220.0
- Positive Control: Kidney Tissue
- Cellular Location: Cell Membrane


*Optimal dilutions should be determined by researchers for the specific applications.

FOR RESEARCH USE ONLY.

AbboMax, Inc 2528 Qume Drive, Suite 8, San Jose, California 95131, USA
1 408-573-1898 (Tel). 1 408-573-1858 (Fax). www.abbomax.com info@abbomax.com



DATA ATTACHMENTS

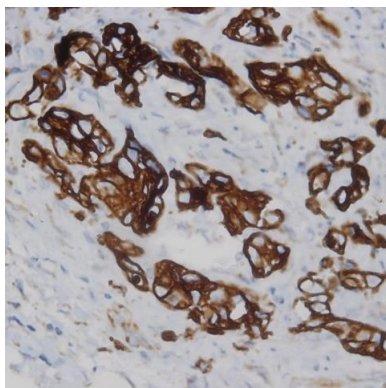

Immunohistochemistry: Human breast carcinoma (FFPE) stained with Mouse anti- Tumor-Associated Glycoprotein (TAG-72) (Cat# 603-880) at 1:200 for 10 min @ RT. Staining of formalin-fixed tissue requires boiling tissue sections in 10 mM Citrate Buffer, pH 6.0 for 10 min followed by cooling at RT for 20 min.

REFERENCES

FOR RESEARCH USE ONLY.

AbboMax, Inc 2528 Qume Drive, Suite 8, San Jose, California 95131, USA
1 408-573-1898 (Tel). 1 408-573-1858 (Fax). www.abbomax.com info@abbomax.com