

Mouse anti CD134 Monoclonal Antibody

Alternative Name(s): TNFRSF4, OX40

Order Information

Description: CD134
Catalogue: 604-150
Lot: See label
Size: 100ug/200ul
Host: Mouse
Clone: L106

• Application: IHC(P), FC

• Reactivity: Hu

ANTIGEN PREPARATION

A recombinant protein of CD134 protein

BACKGROUND

The protein CD134, tumor necrosis factor receptor superfamily, member 4 (TNFRSF4) is a member of the TNF-receptor superfamily. This receptor has been shown to activate NF-kappaB through its interaction with adaptor proteins TRAF2 and TRAF5. Knockout studies in mice suggested that this receptor promotes the expression of apoptosis inhibitors BCL2 and BCL2IL1/BCL2-XL, and thus suppresses apoptosis. The knockout studies also suggested the roles of this receptor in CD4+ T cell response, as well as in T cell-dependent B cell proliferation and differentiation. Studies show that high OX40 (CD134) expression in the tumor immune infiltrate is a favorable prognostic factor of overall survival in non-small cell lung cancer.

PURIFICATION

The mouse IgG is purified by Protein A-Affinity Chromatography according to Isotyping

FORMULATION

This affinity purified antibody is supplied in sterile Phosphatebuffered saline (pH7.2) containing antibody stabilizer

SPECIFICITY

This antibody recognizes human CD134 protein. The other species are not tested.

STORAGE

The antibodies are stable for 24 months from date of receipt when stored at -20oC to -70oC. The antibodies can be stored at 2oC-8oC for three month without detectable loss of activity. Avoid repeated freezing-thawing cycles.

APPLICATIONS/SUGGESTED WORKING DILUTIONS*

• Western Blot: 0.1-1 µg/ml

• ELISA: 0.01-0.1 μg/ml

• Immunoprecipitation: 2-5 µg/ml

• IHC: 2-10 µg/ml

• Flow cytometry: 0.5-5 µg/106 cells


• Molecular Weight: 50.0

• Positive Control: Kidney Tissue

• Cellular Location: Cell Membrane

^{*}Optimal dilutions should be determined by researchers for the specific applications.



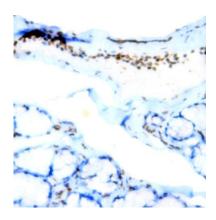

Immunohistochemistry: Human prostate carcinoma (FFPE) stained with Mouse anti-CD134 (Cat# 604-150) at 1:200 for 10 min @ RT. Staining of formalin-fixed tissue requires boiling tissue sections in 10 mM Citrate Buffer, pH 6.0 for 10 min followed by cooling at RT for 20 min.

REFERENCES