

Mouse anti FANCD2 Monoclonal Antibody

Alternative Name(s): FA4; FAD; FACD; FAD2; FA-D2; FANCD

Order Information

Description: FANCD2
Catalogue: 604-830
Lot: See label
Size: 100ug/200ul
Host: Mouse
Clone: FI-17

• Application: IHC(P), WB

• Reactivity: Hu

ANTIGEN PREPARATION

A recombinant protein of human FANCD2

BACKGROUND

Fanconi anemia (FA) is a genetically heterogeneous recessive disorder characterized by cytogenetic instability, hypersensitivity to DNA crosslinking agents, increased chromosomal breakage, and defective DNA repair. The Fanconi anemia complementation group (FANC) currently includes FANCA, FANCB, FANCC, FANCD1 (also called BRCA2), FANCD2, FANCE, FANCF, FANCG, FANCI, FANCJ (also called BRIP1), FANCL, FANCM and FANCN (also called PALB2). FANCD2 protein is monoubiquinated in response to DNA damage, resulting in its localization to nuclear foci with other proteins (BRCA1 AND BRCA2) involved in homology-directed DNA repair.

PURIFICATION

The mouse IaG is purified by Protein A-Affinity Chromatography according to Isotyping

FORMULATION

This affinity purified antibody is supplied in sterile Phosphatebuffered saline (pH7.2) containing antibody stabilizer

SPECIFICITY

This antibody recognizes human FANCD2 protein. The other species are not tested.

STORAGE

The antibodies are stable for 24 months from date of receipt when stored at -20oC to -70oC. The antibodies can be stored at 2oC-8oC for three month without detectable loss of activity. Avoid repeated freezing-thawing cycles.

APPLICATIONS/SUGGESTED WORKING DILUTIONS*

• Western Blot: 0.1-1 μg/ml

• ELISA: 0.01-0.1 μg/ml

• Immunoprecipitation: 2-5 µg/ml

• IHC: 2-10 µg/ml

• Flow cytometry: Not tested


• Molecular Weight: 166.0

Positive Control: Kidney Tissue

• Cellular Location: Cell Membrane

^{*}Optimal dilutions should be determined by researchers for the specific applications.



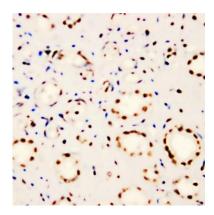

Immunohistochemistry: Human prostate carcinoma (FFPE) stained with Mouse anti- FANCD2 (Cat# 604-830) at 1:200 for 10 min @ RT. Staining of formalinfixed tissue requires boiling tissue sections in 10 mM Citrate Buffer, pH 6.0 for 10 min followed by cooling at RT for 20 min.

REFERENCES