

Mouse anti CD11c Monoclonal Antibody

Alternative Name(s): ITGAX, CD11C; SLEB6

Order Information

• Description: CD11c
• Catalogue: 605-190
• Lot: See label
• Size: 100ug/200ul
• Host: Mouse
• Clone: 3.9

• Application: IHC(P), FC

• Reactivity: Hu

ANTIGEN PREPARATION

A recombinant protein of human CD11c

BACKGROUND

CD11c is a 145-150 kD type I transmembrane glycoprotein. It is also called ITGAX integrin subunit alpha X. Integrins are heterodimeric integral membrane proteins composed of an alpha chain and a beta chain. This I-domain containing alpha integrin combines with the beta 2 chain (ITGB2) to form the integrin lymphocyte function-associated antigen-1 (LFA-1), which is expressed on all leukocytes. CD11c non-covalently associates with integrin β 2 (CD18) and is expressed on monocytes/macrophages, dendritic cells, granulocytes, NK cells, and subsets of T and B cells. CD11c has been reported to play a role in adhesion and CTL killing through its interactions with fibrinogen, CD54, and iC3b.

PURIFICATION

The mouse IgG is purified by Protein A-Affinity Chromatography according to Isotyping

FORMULATION

This affinity purified antibody is supplied in sterile Phosphatebuffered saline (pH7.2) containing antibody stabilizer

SPECIFICITY

This antibody recognizes human CD11c protein. The other species are not tested.

STORAGE

The antibodies are stable for 24 months from date of receipt when stored at -20oC to -70oC. The antibodies can be stored at 2oC-8oC for three month without detectable loss of activity. Avoid repeated freezing-thawing cycles.

APPLICATIONS/SUGGESTED WORKING DILUTIONS*

• Western Blot: 0.1-1 µg/ml

• ELISA: 0.01-0.1 μg/ml

• Immunoprecipitation: 2-5 µg/ml

• IHC: 2-10 µg/ml

• Flow cytometry: 0.5-5 µg/106 cells


Molecular Weight: 130.0

• Positive Control: Kidney Tissue

• Cellular Location: Cell Membrane

^{*}Optimal dilutions should be determined by researchers for the specific applications.



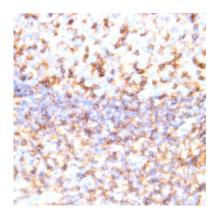

Immunohistochemistry: Human Tonsil (FFPE) stained with Mouse antCD11c (Cat# 605-190) at 1:200 for 10 min @ RT. Staining of formalin-fixed tissue requires boiling tissue sections in 10 mM Citrate Buffer, pH 6.0 for 10 min followed by cooling at RT for 20 min.

REFERENCES