

Mouse anti RRUT-3 Monoclonal Antibody

Alternative Name(s): nan

Order Information

- Description: TIMP-1
- Catalogue: 606-230
- Lot: See label
- Size: 100ug/200ul
- Host: Mouse
- Clone: ZY203/RRUT-3
- Application: IHC(P), WB
- Reactivity: Hu

ANTIGEN PREPARATION

A recombinant protein TIMP-1

BACKGROUND

TIMP-1, TIMP-2, TIMP-3 and TIMP-4 (tissue inhibitor of metalloproteinases 1, 2, 3 and 4) complex with metalloproteinases such as collagenases, gelatin-ases and stromelysins, resulting in irreversible inactivation of the metallopro-teinase. TIMP-1 has been found to be identical to EPA (erythroid-potentiation activity). Parathyroid hormone has been shown to be a regulator of TIMP-2 in osteoblastic cells. TIMP-3 may be involved in regulating trophoblastic invasion of the uterus and remodeling of the extracellular matrix during the folding of epithelia and in the formation, branching and expansion of epithelial tubes. TIMP-4 is most highly expressed in heart, with low levels expressed in liver, brain, lung, thymus and spleen.

PURIFICATION

The mouse IgG is purified by Protein A-Affinity Chromatography according to Isotyping

FORMULATION

This affinity purified antibody is supplied in sterile Phosphatebuffered saline (pH7.2) containing antibody stabilizer

SPECIFICITY

This antibody recognizes human TIMP-1 protein. The other species are not tested.

STORAGE

The antibodies are stable for 24 months from date of receipt when stored at -200C to -700C. The antibodies can be stored at 20C-80C for three month without detectable loss of activity. Avoid repeated freezing-thawing cycles.

APPLICATIONS/SUGGESTED WORKING DILUTIONS*


- Western Blot: 0.1-1 µg/ml
- ELISA: 0.01-0.1 µg/ml
- Immunoprecipitation: 2-5 µg/ml
- IHC: 2-10 µg/ml
- Flow cytometry: Not tested
- Molecular Weight: 71.0
- Positive Control: Kidney Tissue
- Cellular Location: Cell Membrane

*Optimal dilutions should be determined by researchers for the specific applications.

FOR RESEARCH USE ONLY.

AbboMax, Inc 2528 Qume Drive, Suite 8, San Jose, California 95131, USA 1 408-573-1898 (Tel). 1 408-573-1858 (Fax). www.abbomax.com info@abbomax.com



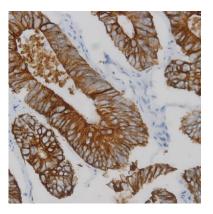

Immunohistochemistry: Human colon carcinoma (FFPE) stained with Mouse anti-TIMP1 (Clone ZY203) (Cat# 606-230) at 1:200 for 10 min @ RT. Staining of formalin-fixed tissue requires boiling tissue sections in 10 mM Citrate Buffer, pH 6.0 for 10 min followed by cooling at RT for 20 min.

REFERENCES