

Mouse anti SSEA-4 Monoclonal Antibody

Alternative Name(s): nan

Order Information

Description: SSEA-4
Catalogue: 606-280
Lot: See label
Size: 100ug/200ul
Host: Mouse
Clone: MC813-70
Application: IHC(P), WB

• Reactivity: Hu

ANTIGEN PREPARATION

A recombinant protein SSEA-4

BACKGROUND

The antibody SSEA-4 (Stage-specific embryonic antigen 4) MC813-70 was raised against human teratocarcinoma cell line. This monoclonal antibody reacts with a glycolipid carbohydrate epitope on SSEA-4 protein. SSEA-4, a cell surface antigen, is a well-known marker of embryonic stem cells, embryonal carcinoma cells and induced pluripotent stem cells. In human cells, these surface markers are down-regulated during the differentiation process. Glycosphingolipids, such as the globo-series GSLs stage-specific embryonic antigen 3 (SSEA-3), SSEA-4, and Globo-H, are known to be associated with various biological processes such as cell recognition, cell adhesion, and signal transduction. SSEA-4 is also frequently overexpressed in cancer cells and cancer stem cells.

PURIFICATION

The mouse IgG is purified by Protein A-Affinity Chromatography according to Isotyping

FORMULATION

This affinity purified antibody is supplied in sterile Phosphatebuffered saline (pH7.2) containing antibody stabilizer

SPECIFICITY

This antibody recognizes human SSEA-4 protein. The other species are not tested.

STORAGE


The antibodies are stable for 24 months from date of receipt when stored at -20oC to -70oC. The antibodies can be stored at 2oC-8oC for three month without detectable loss of activity. Avoid repeated freezing-thawing cycles.

APPLICATIONS/SUGGESTED WORKING DILUTIONS*

- Western Blot: 0.1-1 µg/ml
- ELISA: 0.01-0.1 µg/ml
- Immunoprecipitation: 2-5 µg/ml
- IHC: 2-10 µg/ml
- · Flow cytometry: Not tested
- Molecular Weight: 28.0
- Positive Control: Kidney Tissue
- Cellular Location: Cell Membrane

^{*}Optimal dilutions should be determined by researchers for the specific applications.



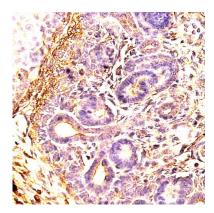

Immunohistochemistry: Mouse embryonic tissue (FFPE) stained with Mouse anti-SSEA4 (Cat# 606-280) at 1:200 for 10 min @ RT. Staining of formalin-fixed tissue requires boiling tissue sections in 10 mM Citrate Buffer, pH 6.0 for 10 min followed by cooling at RT for 20 min

REFERENCES