

Rabbit anti SirT1 Polyclonal antibody

Alternative Name(s): nan

Order Information

Description: SirT1
Catalogue: 630-650
Lot: See label
Size: 100ug/200ul
Host: Rabbit
Clone: nan

• Application: IHC(P), WB • Reactivity: Hu, Ms, Rt,

ANTIGEN PREPARATION

A synthetic peptide NEAISVKQEVTDMNYPSNKS corresponding to human SirT-1. It is identical among human, mice and rat.

BACKGROUND

SIRTs belongs to sirtuin family proteins. The human homologs, which include SIRT1-7, are divided into four main branches: SIRT1-3 are class I, SIRT4 is class II, SIRT5 is class III and SIRT6-7 are class IV. SIRT1 is a member of the sirtuin family of proteins, homologs to the yeast Sir2 protein. Members of the sirtuin family are characterized by a sirtuin core domain and grouped into four classes. SIRT1 is highly expressed in heart, brain and skeletal muscle, with low expression in lung and placenta. SIRT1 regulates the p53-dependent DNA damage response pathway by binding to and deacetylating p53, specifically at Lys 382.

PURIFICATION

The Rabbit IgG is purified by Epitope Affinity Purification

FORMULATION

This affinity purified antibody is supplied in sterile Phosphatebuffered saline (pH7.2) containing antibody stabilizer

SPECIFICITY

This antibody recognizes human SirT1 protein. It cross reacts to human, mice and rat.

STORAGE

The antibodies are stable for 24 months from date of receipt when stored at –20oC to –70oC. The antibodies can be stored at 2oC-8oC for three month without detectable loss of activity. Avoid repeated freezing-thawing cycles.

APPLICATIONS/SUGGESTED WORKING DILUTIONS*

• Western Blot: 0.1-1 μg/ml

• ELISA: 0.01-0.1 μg/ml

• Immunoprecipitation: 2-5 µg/ml

• IHC: 2-10 µg/ml

· Flow cytometry: Not tested


• Molecular Weight: 31.6

• Positive Control: Kidney Tissue

• Cellular Location: Cell Membrane

^{*}Optimal dilutions should be determined by researchers for the specific applications.



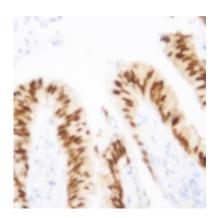

Immunohistochemistry: Human breast carcinoma (FFPE) stained with Rabbit anti-SirT1 antibody (Cat# 630-650) at 1:200 for 10 min @ RT. Staining of formalin-fixed tissue requires boiling tissue sections in 10 mM Citrate Buffer, pH 6.0 for 10 min followed by cooling at RT for 20 min.

REFERENCES